


Main differences between DBLG1 and DBLG2



Differences between DBLG1 and DBLG2

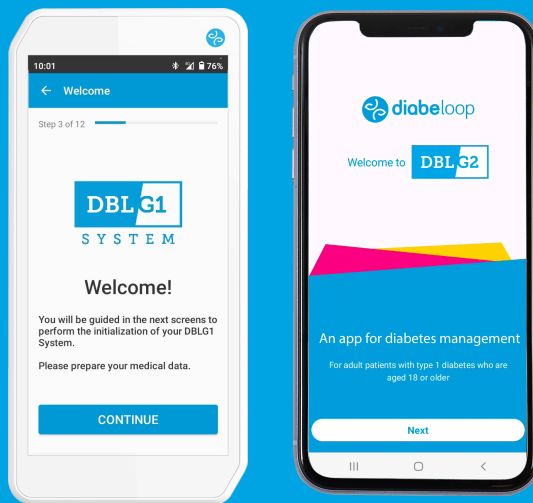

Objective of this memo: to understand the differences between DBLG1 and DBLG2.

DBLG1 is a dedicated handset hosting the Diabeloop algorithm loop mode. DBLG2 is the mobile application version of DBLG1 with some improvements described hereafter.

Starting with the system: downloading and initializing the app



As DBLG2 is a mobile application, you need to download the application and manage updates through the official application store.



Dedicated handset



- **Compatible smartphone** required (list available on Diabeloop's website)



The application is compatible with **Android Operating Systems**.



Your application can be used on **only one smartphone**.

Initialization: addition of extra steps in DBLG2



As DBLG2 is a mobile application, you need to check the sound, grant permissions and activate the app through an activation code.



12 steps



Following additions:

- Check your phone is protected/secure
- Enter an activation code provided by the healthcare professional (or distributor)
- Enter additional information concerning your previous treatment (optional)
- Grant all necessary permissions (e.g., override Do Not Disturb mode, unrestricted battery usage, Bluetooth activation)

! If you change any of these permissions, the application may trigger an alarm or its use may be blocked. Active Internet connection is necessary (via Wi-Fi or mobile data)

Language follows the phone language.

Updates are managed through the official app store.


Some features are managed directly by the phone in DBLG2



Flight mode is managed directly by your phone in DBLG2, like the language. You can use the application in dark mode if you want to.



Flight mode

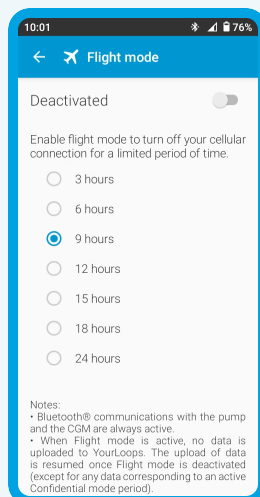

Flight mode

! Some phones deactivate Bluetooth when you activate flight mode. An alarm will be triggered. Remember to reactivate Bluetooth.

Reorganization of the menu and sub-menus



For easier navigation, the app's structure has been reorganized across four pages: Home, Events, System and Profile.



- Rescue carbs
- System status
- History
 - Graphs
 - Events
 - My meals
 - Physical activities
 - Statistics
 - Calibrations
 - Bolus
- Settings
 - Medical
 - Thresholds
 - Loop mode
 - Target glucose level
 - Aggressiveness factors
 - Loop mode OFF (three alerts)
 - Patient
 - Weight & height
 - Typical meals
 - Basal safety profile
 - Total daily dose of insulin
 - Zen mode
 - DBLG1
 - Sound and vibration (alerts/notifications)
 - Credentials
 - PIN Code
 - YourLoops
 - Data sharing
 - Confidential mode
 - Flight mode
 - Languages
 - Preferences (display)
 - Time
 - Updates
 - My profile
- Help: help, support and about



- Home
 - + (declare)
 - Meals
 - Physical activity
 - Rescue carbs
 - Pump bolus
 - Pen bolus
 - Graph
- Events: history of:
 - Declared events (+ edit meals and physical activities when possible),
 - Alarm system events,
 - Boluses (meal and pump/pen)
 - Calibrations
 - Correction boluses (more filters)
- System (= system status)
- Profile
 - Treatment
 - Medical information:
 - Weight, Height
 - Total daily dose of insulin
 - Target and thresholds
 - Aggressiveness
 - Meal size preferences
 - Basal safety profile
 - Insulin type
 - Statistics
 - Notifications
 - Loop mode alerts
 - Loop mode OFF alerts
 - G7 sensor alerts
 - Phone notification settings (Android phones only)
 - Check alarm sound
 - Zen mode
 - Night mode
 - Data sharing
 - Help and contact:
 - Unlock your pump
 - Privacy and security
 - About



Meal management



The rules concerning when you can declare a future meal are different. DBLG2 also gives you the possibility to change the meal type.



Meal

115 mg/dL →
Active Insulin: 1.67 U

Time

-15 min 10:01 +15 min

Quantity Carbs (g)

-5 g 20 +5 g

Small (15g) Medium (20g) Large (25g)

☐ High fat meal

CANCEL CONFIRM

You can declare past and future meals at any time on the same day (from 00:00 to 23:59).



New meal

→ 126 mg/dL Active insulin 1.27 U

Meal details

Start time: Now

Meal type: Breakfast

Meal size

18 g 35 g 53 g Other

☐ High fat

Confirm

You can declare past meals at any time on the same day (from 00:00 to 23:59), but you can only declare future meals one hour ahead.

Physical activity management



The rules concerning when you can declare a physical activity are different.

But both DBLG1 and DBLG2 offer the same list of physical activities to choose from, with the same categorization (aerobic/anaerobic/mixed).



Physical activity

115 mg/dL →

Name: Aerobic (Default)

Start time: Other - Anaerobic

Duration: 15 min

Intensity: Low

Activity: Aerobic (Default)

Other - Anaerobic

Other - Mixed

American football

Badminton

Baseball

Basketball

Combat

Construction work

Cycling

Dancing

You can declare a physical activity at any time on the same day (from 00:00 to 23:59) or for the following day.



New physical activity

→ 125 mg/dL Active insulin 1.27 U

Type

Type: Sport

Time

Start: In 1 hour

Duration: 30 min

Intensity

Low Moderate Intense

Confirm

You can declare a physical in the past or up to 7 days in the future.


Pump bolus management



There are several ways to stop a meal bolus or bolus delivered through the pump control mode.



☰ > System status > Pump section >
MORE > Pump control

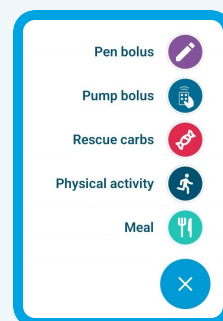

Tap on  > 
or

From the Home screen:

- tap on the "bolus in progress" message. This action takes you to the Pump control screen.
- tap on the bolus or meal pictogram on the graph. This action opens a pop-up window from the bottom of the screen.

From the System screen: go to the pump section and tap More > Pump control or use the shortcut .

From the Events screen: access the bolus directly from the list of current events.

Tap Stop bolus and then Stop bolus again on the pop-up screen. The bolus is updated in the Events list to show the quantity of insulin actually delivered.



External bolus and rescue carbs declarations




Declarations of rescue carbs and external boluses have been grouped with meal and physical activity declarations.



☰ > Rescue carbs

or

☰ > History > Bolus > +

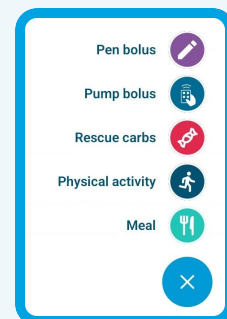

Tap on > for rescue carbs or for an external bolus



Daily insulin consumption on the statistics screen in DBLG2



The statistics screen also displays daily insulin consumption in DBLG2. In DBLG1, this was displayed in ☰ > Settings > Medical > Patient > Total Daily Dose of insulin.

In DBLG2:

- the 1-day tab shows insulin consumption over the last 24 hours.
- the 7-days and 30-days tabs show the average daily insulin over that time period, with daily insulin consumption detailed.

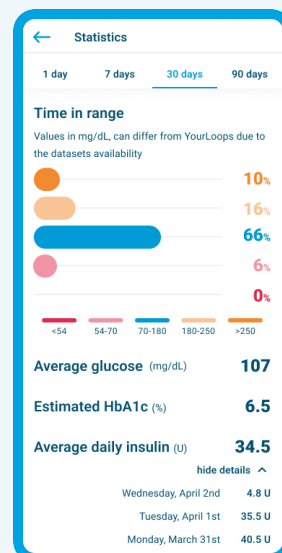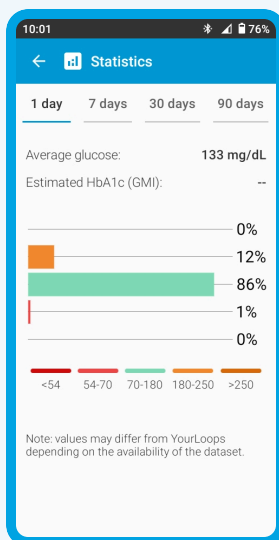


Bluetooth connection request when pairing devices with DBLG2



In DBLG2 you must accept the Bluetooth pairing pop-up when you first pair your devices.

Some pump or sensor statuses have also been reworded in DBLG2.

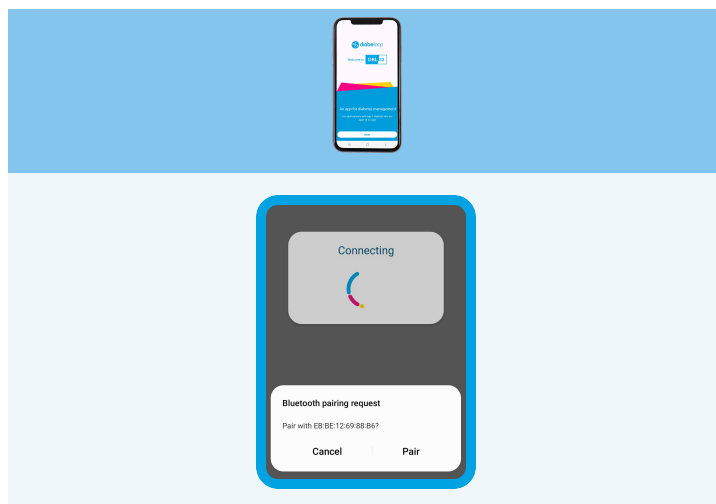


Unlocking the pump with DBLG2



In DBLG2 you can unlock your Kaleido pumps in Profile > Help and contact > Unlock your pump. The steps to follow are described on the screen.

This feature is available since the 1.18 software version of DBLG1.


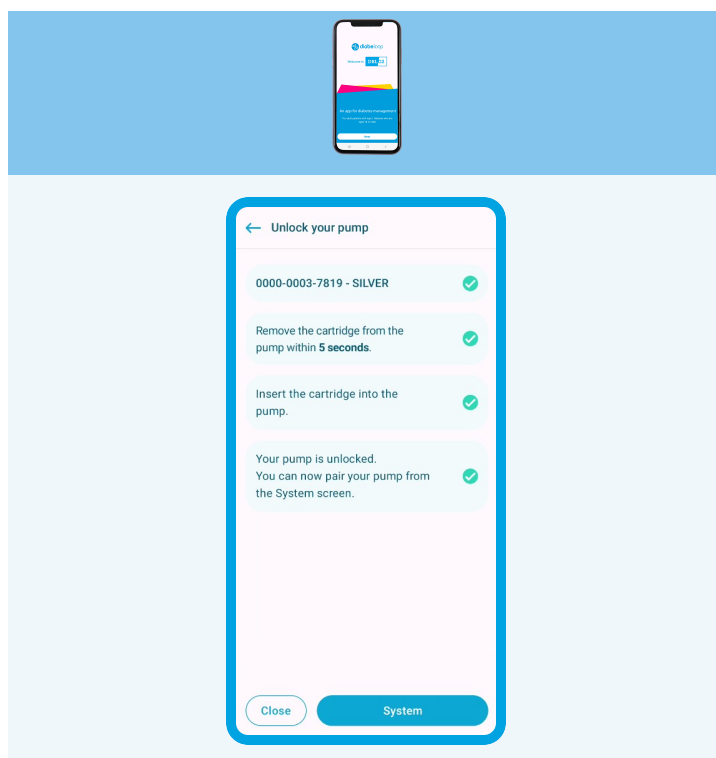


Different rules for manual pump control in DBLG2



The DBLG2 application has a manual control mode, which you can use to:

- deliver or stop a bolus when loop mode is ON or OFF.
- apply (and stop) a temporary basal rate when loop mode is OFF.

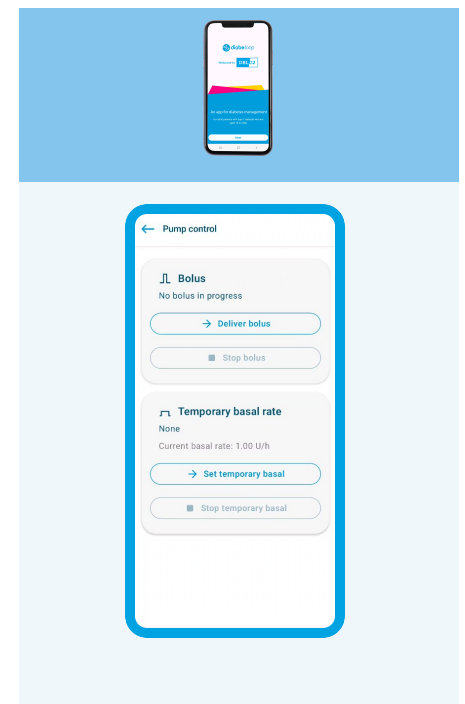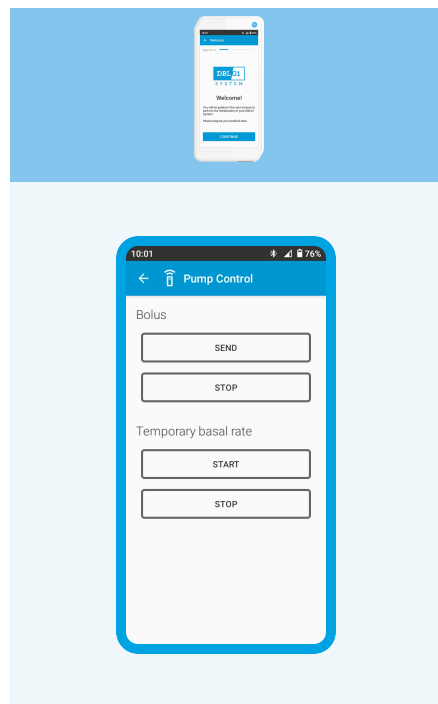


Data transfer in DBLG2 and confidential mode in DBLG1




Your data is automatically pushed to YourLoops in DBLG1 and DBLG2. DBLG1 uses 4G connectivity for data sharing, and you can temporarily stop data sharing via Confidential mode.

DBLG2 uses your mobile connection or Wi-Fi to upload data. You can select Wi-Fi as your preferred method to transfer data. For this, switch on Wi-Fi in Profile > Data sharing > Network preferences.

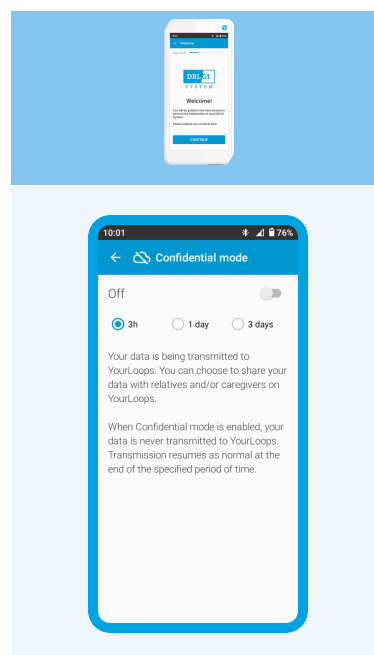

Confidential mode (3 hours, 1 day or 3 days)

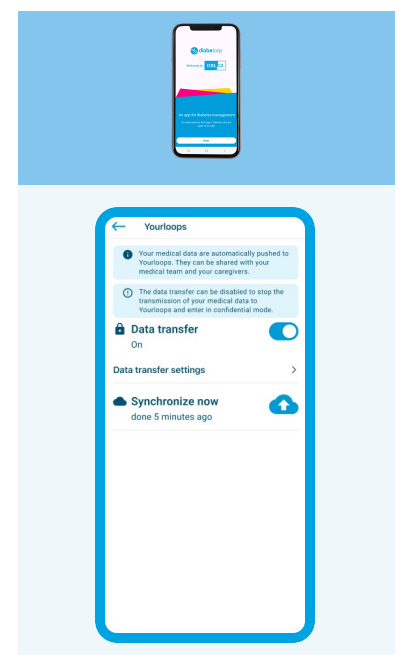

Data transfer to OFF = confidential mode for a specific period of time (3 hours, 12 hours, 1 day or 3 days).



Synchronize now allows you to force synchronization.

Alarm system management through notifications in DBLG2



The sound and vibration of alerts and information messages in DBLG2 follow the phone settings, with the exception of two alerts that can be set to always sound when in loop mode ON.

However, the alarm system works in the same way for DBLG1 and DBLG2.

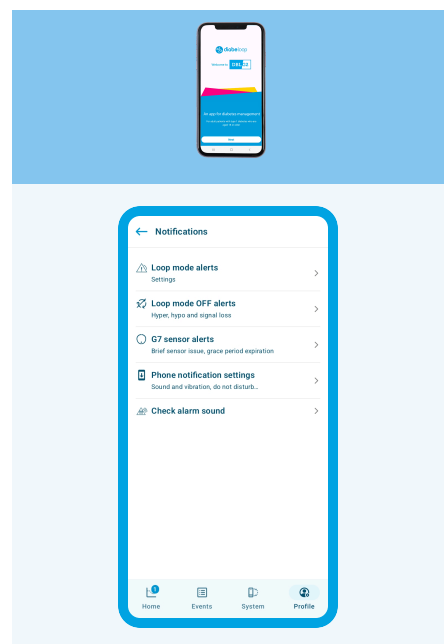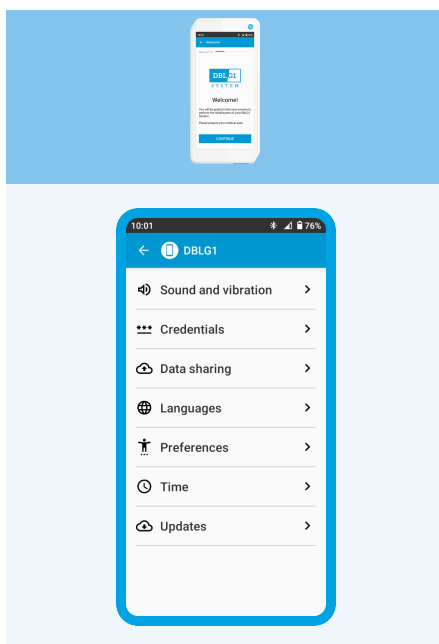


Rescue carbs are recommended (20300) and Rapid rise in glycemia (20102) can be set to always sound, even if your phone is in silent mode.

Night mode for less disturbance at night in DBLG2



Night mode is a new mode in DBLG2 to enable less disturbance at night.



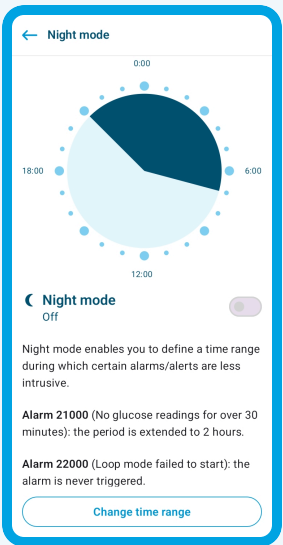

You can define a custom time range for night time. When Night mode is on, these alarms/alerts become less intrusive.

Alarm/alert	In Night mode
21000 No glucose readings for over 30 minutes	Triggered when no readings have been received for over 120 minutes
22000 Loop mode failed to start	Never triggered at night time
20300 Rescue carbs are recommended	The hypoglycemia limit that triggers this alert is lowered by 0.6 mmol/L / 10 mg/dL.

Application management




You can use your activation code on **one smartphone only**. Initializing another smartphone with the same email address and same activation code deactivates the application on the first smartphone.



If you wish to stop using DBLG2, contact your healthcare professional or distributor to have your activation code deactivated. The app will stop automatically afterward.


You will be informed by email three days before the deactivation, and a temporary message will appear in your app showing the exact end date. When that date is reached, you will be prompted to unpair and unlock your pump(s). Once these steps have been completed, you can safely uninstall the app from your phone. Your personal and medical data will be deleted from your phone.

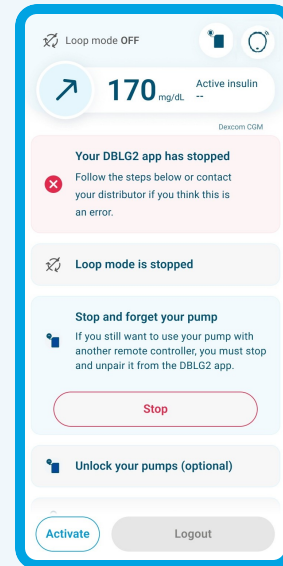

If you want to use your DBLG2 application on a new phone, go to **Profile > My account > Manage account** and log out from your account.

Then install the application on the new phone and initialize it using the same activation code, your YourLoops credentials and the settings visible on YourLoops, or with the help of your healthcare professional.




We advise you to not log out if you don't want to change your phone, as this will delete all the application data and you will have to start all over again from the initialization phase. Moreover, you won't be able to reconnect to YourLoops if you delete your account.

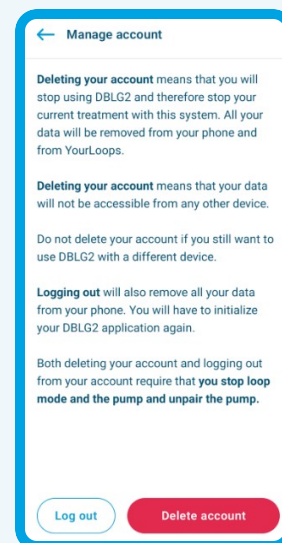

Help and contact



In DBLG2 you can access an electronic version of the user guide and a webhelp in the Help section.

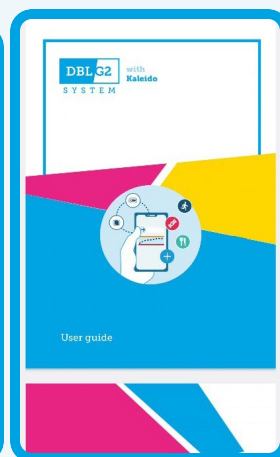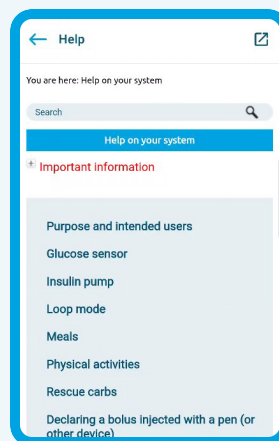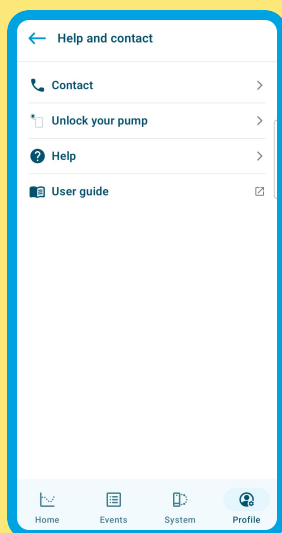

History

DA Number	Version	Description	Effective date	Author
DA 00178	0	Creation	2025-04-28	Auriane VICENTE
DA 00178	1	Remove of DBLG2 algorithm equivalent to software version 1.16 of DBLG1 (now equivalent to 1.17 software version). Small wording correction.	See Qualios approval record	Auriane VICENTE
DA 00268	2	Updated to align with Version android 1.2.0(G7 commercial)	See Qualios approval record	Cynthia BACSO



DBLG2-ED-ME-001-en release 2

Main differences between DBLG1 and DBLG2

Nature of changes :

Updated to align with Version android 1.2.0(G7 commercial)

Attached files :

DBLG2-ED-ME-001-en-Rev2-Main differences between DBLG1 and DBLG2.pdf

Status : Active

Submitted on : 13/02/2026

Approved on : 13/02/2026

Activated on : 13/02/2026

Next review on : 13/02/2028

Authors :

BACSO Cynthia - Education Specialist

Editors :

DIEUX Lorine - QA Specialist

Signatory :

1. Author - BACSO Cynthia - Education Specialist digitally signed on 13/02/2026 11:12
2. Reviewer - CARON Martine - Technical Documentation Specialist digitally signed on 13/02/2026 14:18
3. Approver - DIEUX Lorine - QA Specialist digitally signed on 13/02/2026 14:20

Distribution groups :

Tous les utilisateurs actifs - E-Mail

Departments in distribution :

Education - E-Mail

Users in distribution :

Bot Qualios Slack - E-Mail

Document history :

DBLG2-ED-ME-001-en release 2 : Main differences between DBLG1 and DBLG2

Approved on 13/02/2026

Activated on 13/02/2026

Nature of changes : Updated to align with Version android 1.2.0(G7 commercial)

DBLG2-ED-ME-001-en release 1 : Main differences between DBLG1 and DBLG2

Approved on 05/06/2025

Activated on 05/06/2025

Nature of changes : See doc history

Archived on 13/02/2026

Consulted on : 16/02/2026 09:05